




BOOKWALTER™ RETRACTOR KIT II

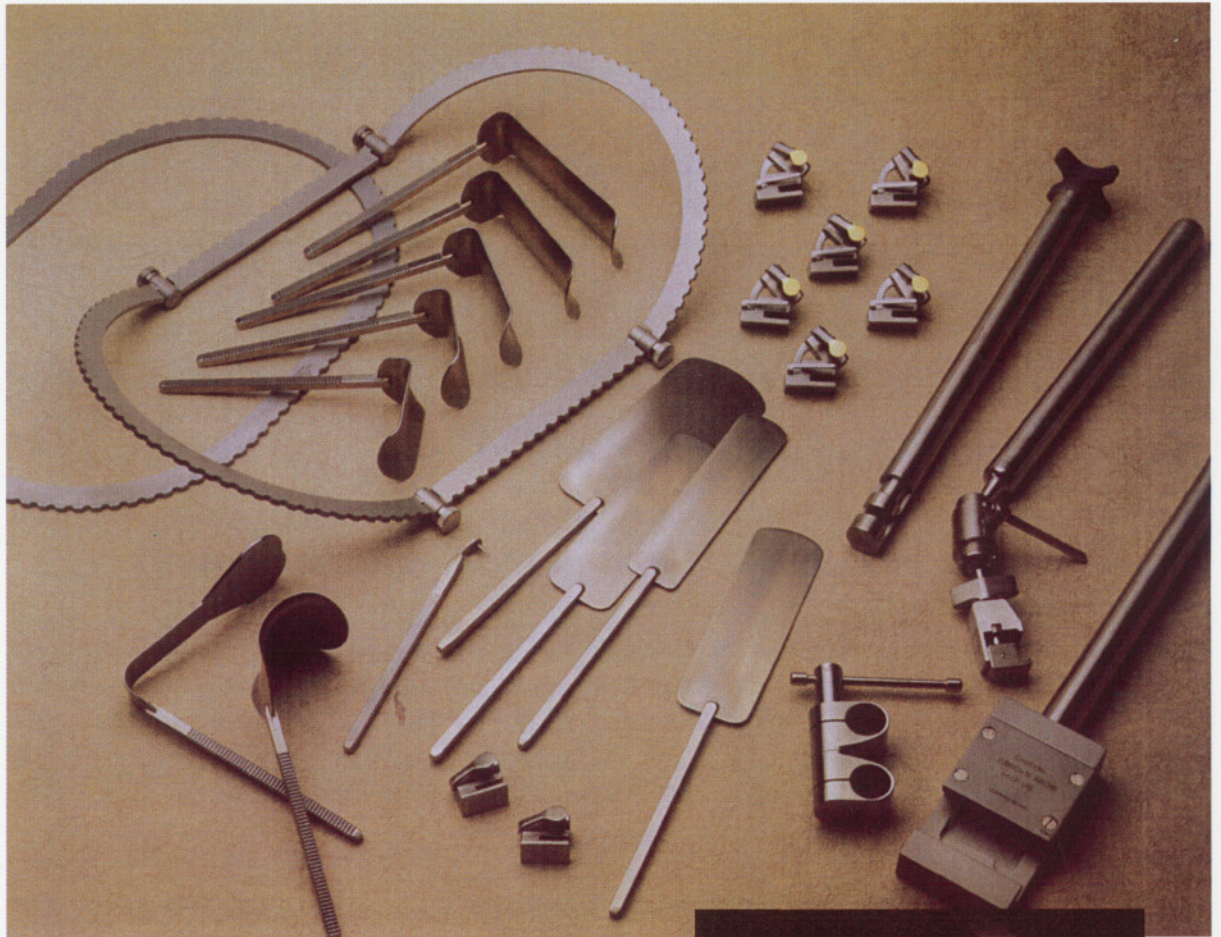


The BOOKWALTER Fiberoptic Light fits securely to any BOOKWALTER ring. With its memory sleeve construction, it can be contoured to afford the surgeon direct lighting during the surgical procedure.

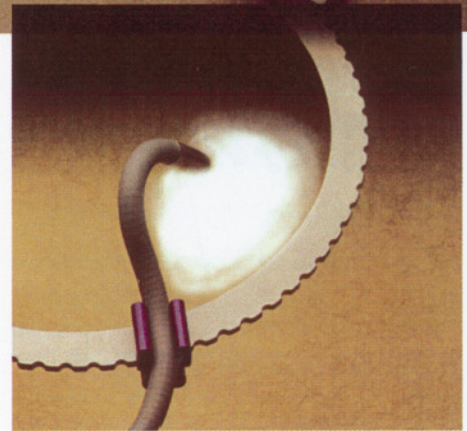


Codman

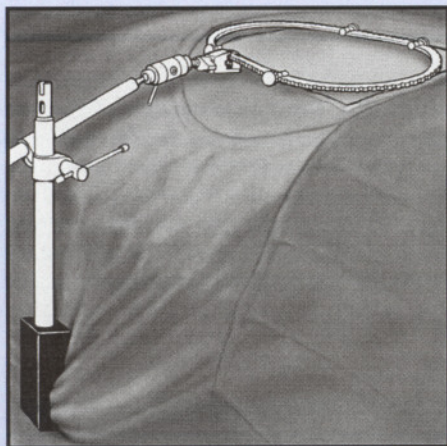


The appropriate ring (medium oval and segmented are included in kit) is selected, attached to the horizontal bar (or horizontal flex bar) and centered over the operative site, not necessarily over the incision.

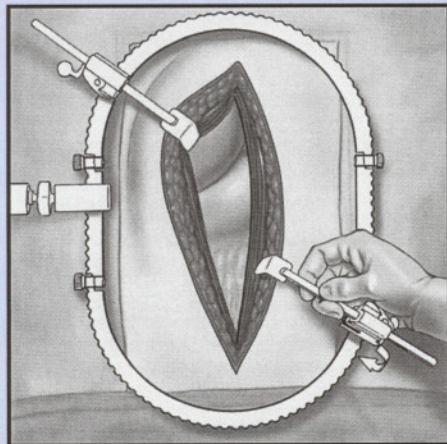

Select the appropriate retractor blades, and attach them to the ring via the quick release ratchet mechanisms.

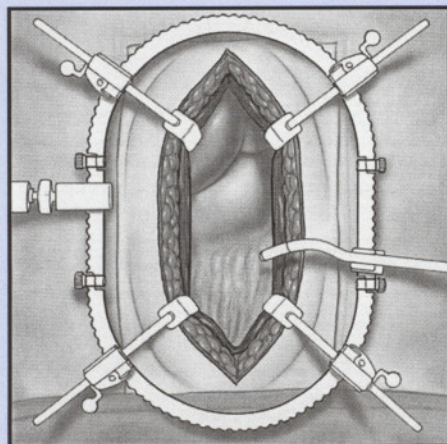

The blades are positioned to provide optimal exposure. Readjustment, repositioning or replacement is fast and simple.

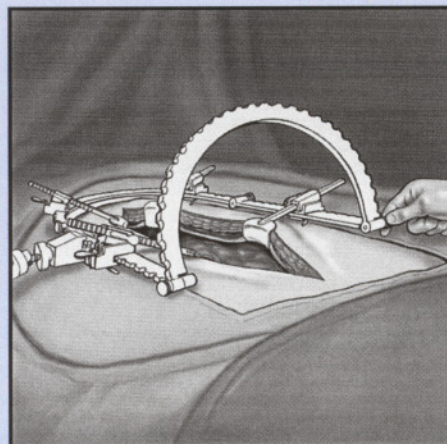

The BOOKWALTER Retractor provides simple, multi-directional exposure for a broad range of surgical procedures.


BOOKWALTER is a trademark of John Bookwalter, M.D.

©2006 Codman & Shurtleff, Inc. All rights reserved.

BOOKWALTER Retractor Kit II

The BOOKWALTER Kit II can be used effectively for many procedures. **Single post fixation**, choice of **ring configuration**, **retractor blade options**, **quick release ratchet mechanism**, and choice of **horizontal bar** or **horizontal flex bar** combine to afford optimal exposure important in many procedures, including cholecystectomy, retroperitoneal resections, nephrectomy, abdominal aortic aneurysms, iliac femoral grafts, abdominal hysterectomy, appendectomy, hernia repair and transplant procedures.

ORDERING INFORMATION

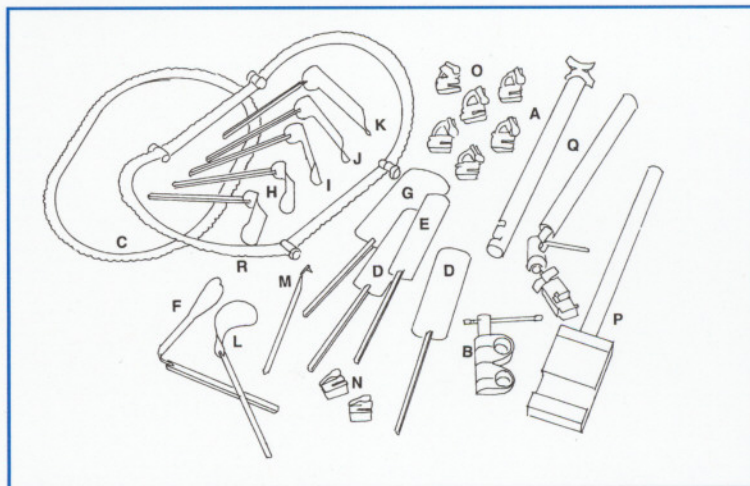

50-4542

BOOKWALTER Retractor Kit II

Kit Includes:

- | | |
|------------|--|
| A. 50-4552 | Horizontal Bar |
| B. 50-4554 | Post Coupling |
| C. 50-4557 | Medium Oval Ring, 10 1/2 X 14 1/2" I.D. (26.7cm X 36.8cm) |
| D. 50-4561 | Malleable Retractor Blade, 2 X 4" (51mm X 10.1cm) |
| E. 50-4562 | Malleable Retractor Blade, 1 1/2 X 6" (38mm X 15.2cm) |
| F. 50-4566 | Harrington Retractor Blade |
| G. 50-4567 | Malleable Retractor Blade, 3 X 6" (76mm X 15.2cm) |
| H. 50-4568 | Kelly Retractor Blade (2 ea.), 2X 3" (51mm X 76mm) |
| I. 50-4569 | Kelly Retractor Blade, 2 X 4" (51mm X 10.1cm) |
| J. 50-4570 | Kelly Retractor Blade, 2 X 5" (51mm X 12.7cm) |
| K. 50-4571 | Kelly Retractor Blade, 2 X 6" (51mm X 15.2cm) |
| L. 50-4576 | Balfour Retractor Blade |
| M. 50-4577 | Gelpi Retractor Blade |
| N. 50-4579 | Ratchet Mechanism (2 ea.) |
| O. 50-4580 | Tilt Ratchet Mechanism (6 ea.) |
| P. 50-4581 | Table Post |
| Q. 50-4582 | Horizontal Flex Bar |
| R. 50-4583 | Medium Segmented Ring, 10 1/2 X 19 1/2" I.D. (26.6cm x 49.5cm) |

Also Available:

- | | |
|---------|---|
| 50-4627 | Fiberoptic Light Set (Requires Fiberoptic Cable with standard 5mm bundle) |
|---------|---|

Set Includes:

- | | |
|---------|--------------------------------------|
| 50-4625 | Fiberoptic Clip |
| 50-4626 | Fiberoptic Light (approx. 12", 30cm) |

Codman
a Johnson & Johnson company

For more information, contact your Codman sales representative

325 PARAMOUNT DRIVE • RAYNHAM, MA 02767 • (508) 880-8100
FOR PRODUCT INFORMATION, CALL 1-800-225-0460

RET-09-000 03/06 QP/QP